




CIRCA
SCIENTIFIC

CrossWise™ RF



Transseptal Access System

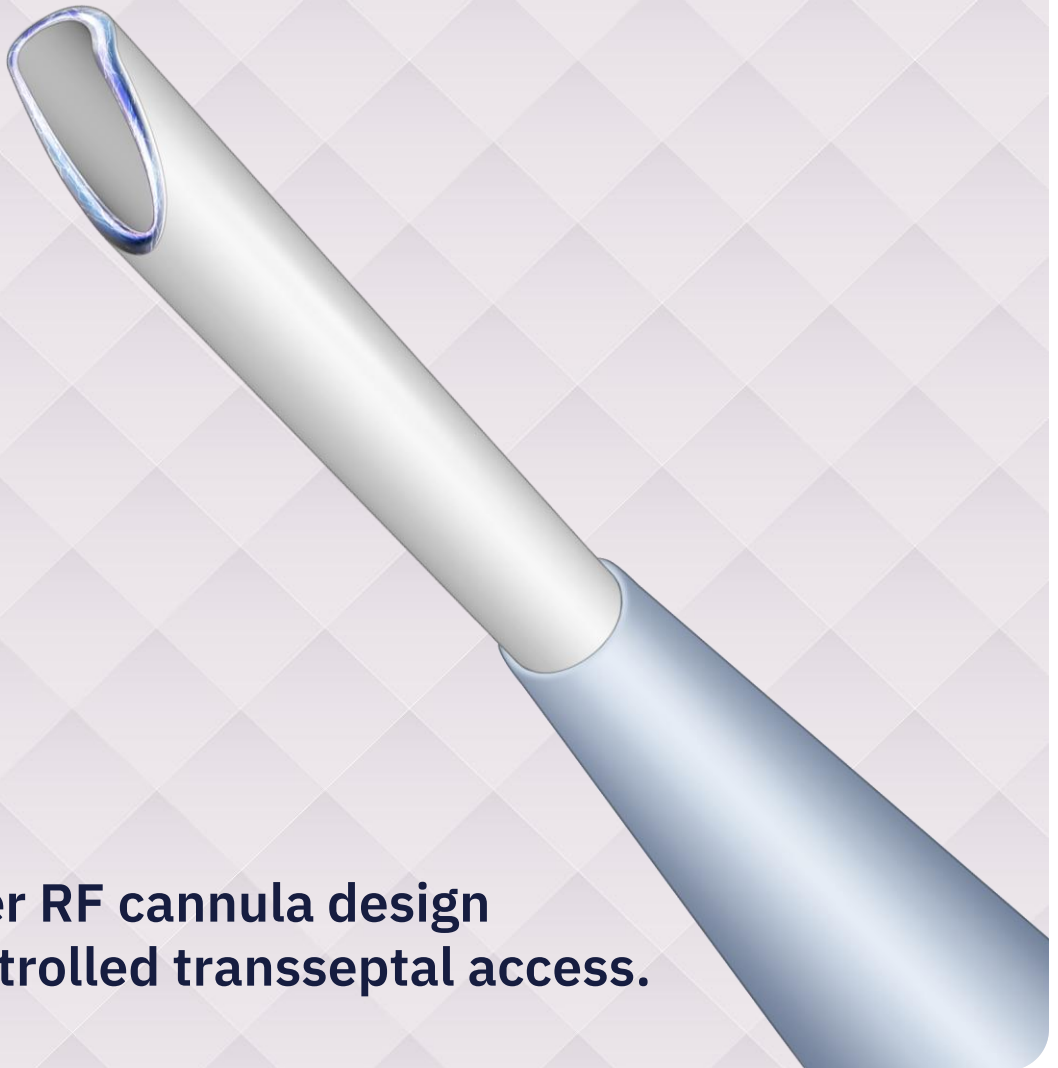


**Smarter RF cannula design
for controlled transseptal access.**

Controlled Crossing

The CrossWise™ RF Transseptal Access System is engineered for precise control. The novel cannula design provides a stable crossing even in the most challenging anatomies.

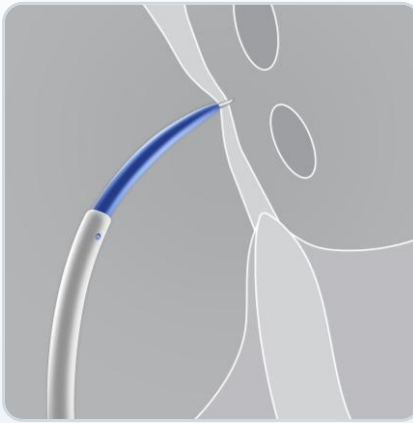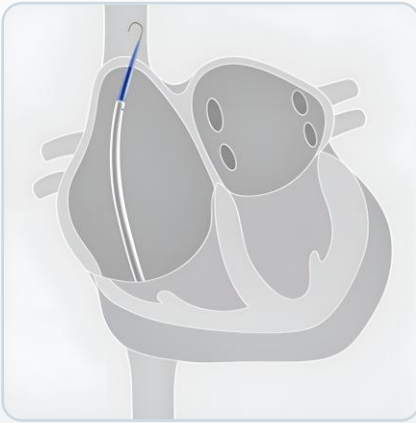



Patented atraumatic beveled tip

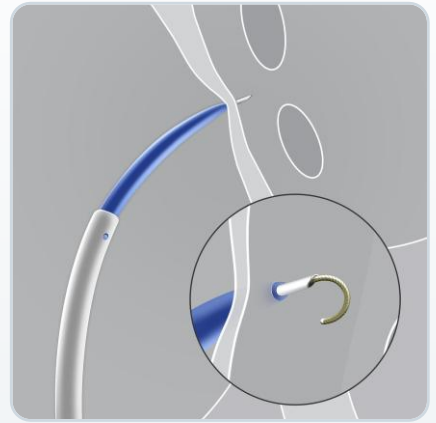

Quick clean puncture

Visual Safety Indicator

Sleeve provides clear visual confirmation when the cannula tip is protected within the dilator, helping to ensure controlled and deliberate advancement.

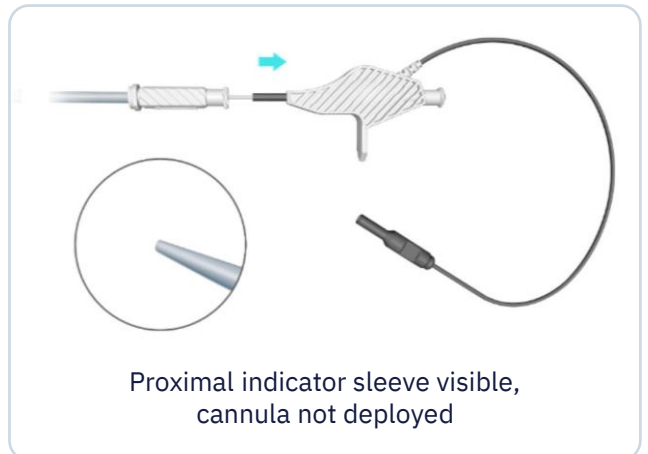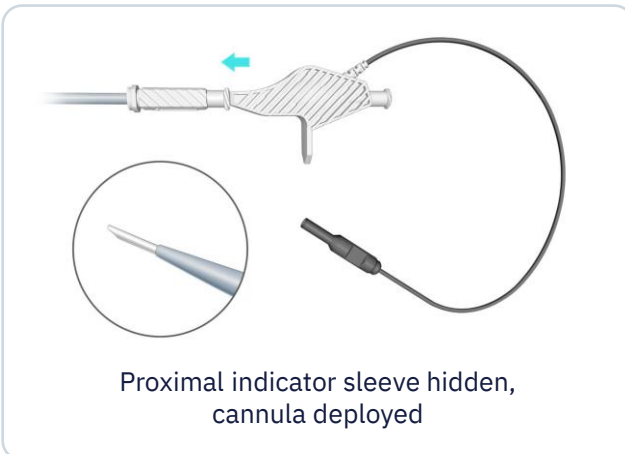

Confidence

- Cannula designed for clean consistent puncture
- The crescent-shaped tissue defect is created without the risk of coring¹
- Successful puncture achieved in as little as 290ms¹



Choice

Compatible with commonly used steerable or fixed sheaths with no proprietary capital equipment required.* Use with the tools and techniques of your preferred workflow.^

Connector to generator system



Indicator Sleeve

Shapable Dilator

0.032" J Tip
Guidewire

Dielectrically
coated cannula

Zero Exchanges

Cannula and dilator
are removed

Wire stays in place
for therapy



Our Cardiocentric™ Mission

To advance procedural techniques and deliver innovative cardiac solutions that help physicians reduce complications and improve procedural success.

Ordering Information

| Product Code | Description |
|-----------------|--|
| CW-1085A | CrossWise RF Transseptal Access System – 8.5Fr, 97CM, C1, Agilis† |
| CW-1085S | CrossWise RF Transseptal Access System – 8.5Fr, 71CM, C1, Swartz† |
| CW-1012C | CrossWise RF Transseptal Access System – 12Fr, 88cm, C1, FlexCath† |
| CW-1012W | CrossWise RF Transseptal Access System – 12Fr, 88cm, C1, Watchman† |
| CW-1013F | CrossWise RF Transseptal Access System – 13Fr, 97CM, C1, Faradrive† |
| CW-1085V | CrossWise RF Transseptal Access System – 8.5Fr, 97CM, C1, Vizigo† |
| CW-1001 | CrossWise Multi-Use RF Adapter Cable Connects CrossWise to commercially available electrosurgical pencil |



CIRCA Scientific, Inc.
 Corporate Office
 14 Inverness Drive East
 Suite H-136
 Englewood, CO 80112



info@circascientific.com



1.303.951.8767



www.circascientific.com

Indications for Use

The CrossWise™ RF Transseptal Cannula and accessories are used to create an atrial septal defect in the heart. Secondary indications include infusing solutions including heparinized saline and mixtures of 50% contrast media and 50% saline.

RX Only Caution

Federal (U.S.A.) law restricts the use of this device for sale by or on the order of a physician.

Patents: www.circascientific.com

© 2026 CIRCA Scientific, Inc.

All rights reserved.

CrossWise is a trademark of CIRCA Scientific, Inc.† The third-party trademarks used herein are trademarks of their respective owners.

CW-ART 100 Rev. 04

¹ Ibrahim R, Hasselwander C, Lam S, Lloyd M. Validation of a core-less transseptal system using a specialized energized cannula. Poster presented at: Heart Rhythm Society (HRS) Scientific Sessions 2025; April 30–May 3, 2025; San Diego, CA.

*Generator is not included with CrossWise RF Transseptal Access System.

^Refer to IFU for full instructions.